



Hillrom™



Welch Allyn®
ABPM 7100 Ambulatory
Blood Pressure Monitor

PRODUCT SPECIFICATIONS + PARTS AND ACCESSORIES

Devices	Description
ABPM-7100	ABPM 7100 Recorder
ABPM-7100S	ABPM 7100 Recorder including CardioPerfect WorkStation Software
ABPM-7100CBP	ABPM 7100 Recorder with Central Blood Pressure Option and HMS Software

Accessories	Description
7100-21	ABPM 7100 Pouch & Shoulder Belt
7100-24	ABPM 7100 USB Interface Cable
REUSE-10-ABPM	ABPM Cuff Small Adult: 20 – 24 cm
REUSE-11-ABPM	ABPM Cuff Adult: 24 – 32 cm
REUSE-11L-ABPM	ABPM Cuff Adult Plus: 32 – 38 cm
REUSE-12-ABPM	ABPM Cuff Extra Large: 38 – 55 cm
REUSE-US-ABPM	Cuff Set (10, 11, 11L, 12 sizes)

For more information, please contact your local distributor or Hillrom sales representative at 1-800-445-3730 or learn more at www.hillrom.com.

hillrom.com

4341 State Street Rd., Skaneateles Falls, NY 13153

¹ Blood Pressure & Pulse Wave Measurement combined in one procedure Re-classification of risk patients: I.E.M. GmbH, 52222 Stolberg, Germany.
Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.



Hillrom™



Welch Allyn®
ABPM 7100 Ambulatory
Blood Pressure Monitor

Comprehensive blood pressure monitoring from the comfort of home

Getting an accurate picture of your patient's blood pressure is crucial to making the right call for their care. The Welch Allyn® ABPM 7100 Ambulatory Blood Pressure Monitor is an easy-to-use, 24-hour BP monitor that's designed to help you do just that.

Monitor your patients outside the office to help avoid the effects of white coat hypertension and take accurate sleep readings so you can better tailor drug therapy regimens to meet patient needs. Our powerful software can help you enhance patient blood pressure management while making at-home monitoring more efficient for your practice.




Powerful software helps make data management easier

The CardioPerfect WorkStation software suite manages data from the ABPM 7100 and other Welch Allyn cardiology devices in one central system, allowing you to store all patient information in the same record with paperless efficiency. Retrieve and analyze information quickly with a simple interface, then send patient information to your EMR — so you can save time and focus on delivering quality patient care.

EASILY CREATE, PRINT AND CUSTOMIZE DETAILED REPORTS:

- 1 Statistical summaries for systolic, pulse, mean arterial pressure and heart rate (including minimums, maximums, average and standard deviation)
- 2 Graphical representation of a patient's BP profile (including pie charts, histograms, correlation charts and day vs night reading breakdowns)
- 3 Editable, tabular report for detailed analysis
- 4 Trending tool for long-term evaluation and help with drug management

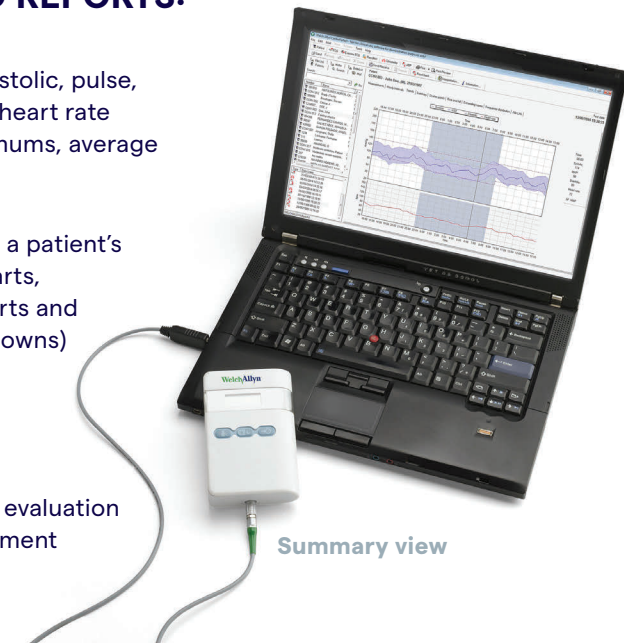

Experience comprehensive blood pressure analysis

- Unlimited protocols allow you to program measurement periods and inflation frequencies to meet any study needs
- Auto Feedback Logic reduces over-inflation and provides improved patient comfort and acceptance
- The lightweight, compact and quiet design helps improve patient compliance
- Welch Allyn CardioPerfect® WorkStation software helps you analyze readings to diagnose potential complications
- Meets or exceeds the latest AAMI, EHS Standards and has a BHS A/A Rating



The benefits of measuring central blood pressure¹

Studies show that measuring aortic central blood pressure provides a greater indication of potential heart attack and stroke than a blood pressure taken on the upper arm alone. The ABPM 7100 has an upgradable Central Blood Pressure option that uses proprietary technology to help you more accurately detect cardiovascular risk and subclinical organ damage.

