



Hillrom™

Welch Allyn® Green Series™
Medical Exam Lighting



BRIGHT LIGHT WHEREVER YOU NEED IT, FOR LESS

Our Green Series Medical Exam Lights give off a brighter, whiter LED light and consume less energy than traditional halogen lights. See true tissue color in a well-defined lit field while also reducing your cost of ownership.



BRIGHTER

- View true tissue color with 5,500° K
- Get a focused, uniform light with crisp edges and no dark or hot spots



GREENER

- Lower costs with long-lasting LEDs that do not have to be replaced
- Consume less energy

Welch Allyn® Green Series™ Medical Exam Lights



GS Exam
Light IV

GS EXAM LIGHT IV

- Ideal for OB/GYN offices or specialties that require focused light in specific areas
- One LED with 3X the typical light output of halogen
- Get pinpoint accuracy with a compact, unobtrusive head and a long, flexible light pipe



GS 300
General
Exam Light

GS 300 GENERAL EXAM LIGHT

- Ideal for primary care facilities, emergency rooms, outpatient facilities and dermatology
- One LED gives a perfect combination of intense light output and cool operation, all delivered in a compact design
- Affordable price and rugged performance

MORE PATIENT-FRIENDLY

- Reduce risk of cross-contamination with touchless controls and Welch Allyn® Disposable Sheaths
- Improve patient comfort and safety with minimal heat dissipation



Welch Allyn® Green Series™ Procedure Lights



GS 600 Minor
Procedure Light

GS 600 MINOR PROCEDURE LIGHT

- Ideal for minor procedures and exams in any healthcare environment
- Three LEDs and a wider head give exceptional light in demanding viewing conditions
- Adjust for optimal output with touchless intensity controls



GS 900
Procedure Light

GS 900 PROCEDURE LIGHT

- Ideal for procedures in your hospital, ambulatory care center, outpatient facility or physician's office
- Six LEDs with intensity controls and outstanding durability, maneuverability and spot quality
- Maximize floor space and get direct light wherever you need it with ceiling mount option

Welch Allyn® Green Series™ Medical Headlights



GS PORTABLE HEADLIGHTS

- Ideal for general examination, minor suturing and foreign body removal
- Enjoy maximum comfort with a weight of only 4 ounces and choice of rigid or terrycloth headband
- Rechargeable batteries within headband give 60 minutes of continuous runtime from full charge
- LED provides up to 10,000 hours of life
- Coaxial luminaire gives well-defined, shadow-free illumination



GS PROCEDURE HEADLIGHT

- Ideal for any procedure or application that requires a small, bright light, such as foreign body removal, OB/ GYN and dermatology
- Enjoy maximum comfort with a weight of only 6 ounces
- Use with direct power supply or a rechargeable battery that gives 60 minutes of continuous runtime from full charge
- LED provides up to 50,000 hours of life
- Coaxial luminaire gives well-defined, shadow-free illumination

When you need targeted, well-defined illumination, look to our high quality, long-lasting and reliable Green Series Procedure Headlights to improve staff satisfaction and enhance patient care.



Specifications



EXAM LIGHTS

	GS Exam Light IV	GS 300 General Exam Light
Number of LEDs	1 White LED	1 White LED
Color Temperature	5,500° K	5,500° K
Light Intensity at Typical Working Distance	20,000 lux (1,858 footcandles)	12" distance: 30,000 lux (2,787 footcandles) 20" distance: 15,000 lux (1,393 footcandles)
LED Life	50,000 hours	50,000 hours
Working Distance	16 in (40 cm)	12 in (30.5 cm) to 20 in (50 cm)
Spot Size at Working Distance	4 in (10 cm) to 10 in (25 cm)	5 in (13 cm) to 6 in (15 cm)
Electrical	Input: 100-240 V~ 50 Hz-60 Hz, 0.30 A max	Input: 100-240 V~ 50 Hz-60 Hz, 0.35 A max
Power Cord	#18 AWG, 3-wire grounded cord set, 8' (2.5 m) long	#18 AWG, 3-wire grounded cord set, 8' (2.5 m) long
Weight	Total Table/Wall Mount: 5.30 lbs (2.40 kg) Total Mobile Stand: 12.40 lbs (5.63 kg)	Total Table/Wall Mount: 5.65 lbs (2.56 kg) Total Mobile Stand: 12.75 lbs (5.78 kg)
Mounting Options	Table/Wall Mount Mobile Stand	Table/Wall Mount Mobile Stand
Accessory Options	Disposable Sheath #52640	NA
Warranty	5-year limited	5-year limited



PROCEDURE LIGHTS


HEADLIGHTS

GS 600 Minor Procedure Light	GS 900 Procedure Light	GS Portable Headlights	GS Procedure Headlight
3 White LEDs	6 White LEDs	LED	5,700° K
5,500° K	5,500° K	5,300° K	5,700° K
High Intensity Setting: 12" distance: 75,000 lux (6,967 footcandles) 20" distance: 40,000 lux (3,716 footcandles) Low Intensity Setting: 12" distance: 55,000 lux (5,109 footcandles) 20" distance: 30,000 lux (2,787 footcandles)	High Intensity Setting: 24" distance: 100,000 lux (9,290 footcandles) 40" distance: 55,000 lux (5,109 footcandles) Low Intensity Setting: 24" distance: 75,000 lux (6,967 footcandles) 40" distance: 40,000 lux (3,716 footcandles)	2,000 lux, 40 lumens	9,000 lux, 120 lumens
50,000 hours	50,000 hours	10,000 hours	18 in (46 cm)
12 in (30.5 cm) to 20 in (50 cm)	24 in (60.9 cm) to 40 in (100 cm)	16 in (40 cm)	18 in (46 cm)
6 in (15 cm) to 7 in (18 cm)	9 in (23 cm) to 10 in (25 cm)	4 in (10 cm) to 8 in (20 cm)	1 in (2.5 cm) to 6 in (15 cm)
Input: 100-240 V~ 50 Hz-60 Hz, 0.65 A max	Input: 100-240 V~ 50 Hz-60 Hz, 1.00 A max	100-240 V~ 50 Hz-60 Hz, 0.2 A max	100-240 V~ 50 Hz-60 Hz, 0.2 A max
#18 AWG, 3-wire grounded cord set, 8' (2.5 m) long	#18 AWG, 3-wire grounded cord set, 8' (2.5 m) long		
Total Table/Wall Mount: 6.45 lbs (2.93 kg) Total Mobile Stand: 17.35 lbs (7.87 kg)	Total Wall Mount: 17.90 lbs (8.12 kg) Total Mobile Stand: 44.85 lbs (20.34 kg) Total Ceiling Mount: 4.50 lbs (11.11 kg)	4 ounces (113 g)	6 ounces (170 g)
Table/Wall Mount Mobile Stand	Wall Mount Mobile Stand Ceiling Mount	Head Mount-Terrycloth or Rigid Headband	Head Mount-Rigid Headband
Disposable Sheath #52630	Disposable Sheath #52630	Soft Case #05232-U	Carrying Case #49099, Battery Pack with Charger #75260
5-year limited	5-year limited	1-year limited	1-year limited



Hillrom™

CHOOSE FROM A VARIETY OF MOUNTING OPTIONS

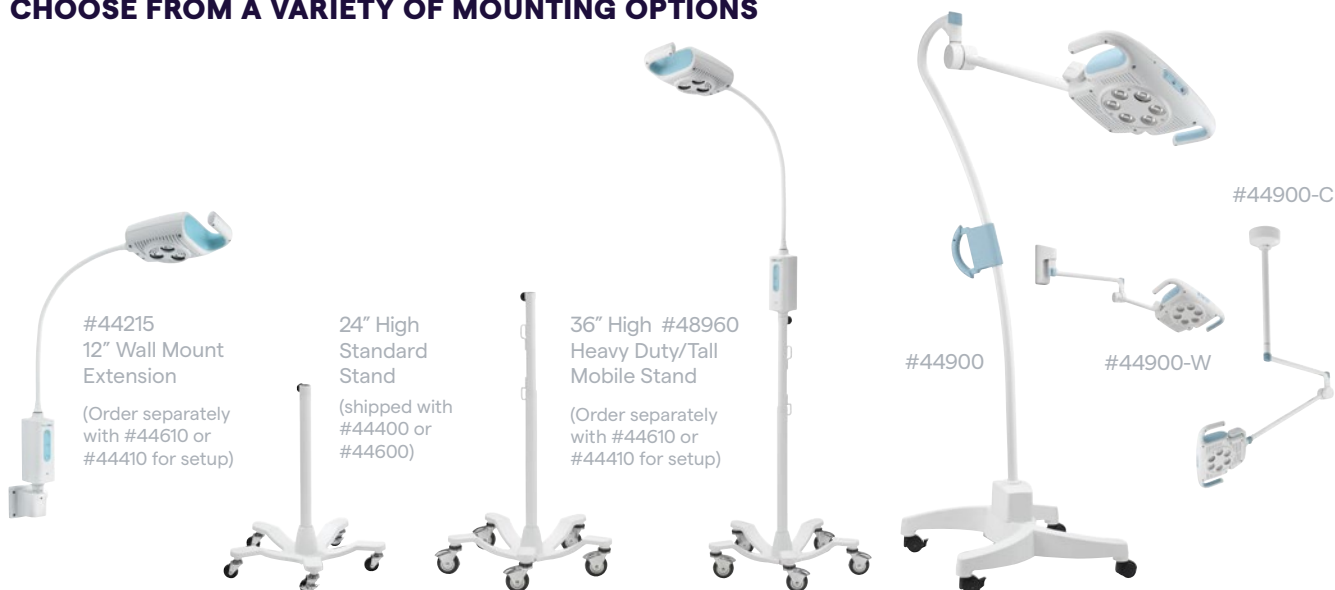

PRODUCT AND ACCESSORY ORDERING INFORMATION

GS EXAM LIGHT IV

48810: GS Exam Light IV with Table/Wall Mount
48950: Mobile Stand for GS Exam Light IV
52640: Disposable Sheath for GS Exam Light IV
(25/box, 5 boxes/case)

GS 300 GENERAL EXAM LIGHT

44400: GS 300 General Exam Light with Mobile Stand
44410: GS 300 General Exam Light with Table/Wall Mount
44215: 12" Wall Mount Extension
48956: Exam Light IV Midmark Bracket Adapter for Midmark Tables (Used to mount Exam Light IV to Midmark Spline Mount Bracket)
48960: Heavy Duty/Tall Mobile Stand

GS 600 MINOR PROCEDURE LIGHT

44600: GS 600 Minor Procedure Light with Mobile Stand
44610: GS 600 Minor Procedure Light with Table/Wall Mount
44215: 12" Wall Mount Extension
48960: Heavy Duty/Tall Mobile Stand
52630: Disposable Sheath for GS 600/GS 900
(50/box, 5 boxes/case)

GS 900 PROCEDURE LIGHT

44900: GS 900 Procedure Light with Mobile Stand
44900-C: GS 900 Procedure Light with Ceiling Mount
44900-W: GS 900 Procedure Light with Wall Mount
52630: Disposable Sheath for GS 600/GS 900
(50/box, 5 boxes/case)

GS PORTABLE HEADLIGHT

46070R: GS Portable Headlight with Rigid Headband and Direct Power Supply/Charger
46070: GS Portable Headlight with Terrycloth Headband and Direct Power Supply/Charger

GS PROCEDURE LIGHT

49020: GS Procedure Headlight with Direct Power Supply
75260: GS Procedure Headlight Battery Pack with Charger

For more information, please contact your local distributor or Hillrom sales representative at 1-800-445-3730.

hillrom.com

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

© 2021 Welch Allyn, Inc. ALL RIGHTS RESERVED. APR129104 Rev 1 01-JUL-2021 ENG – US